## ETB 100 SONIC TOOTHBRUSH USER MANUAL

### PRODUCT INTRODUCTION

The sonic toothbrush is a high-frequency vibration type electric toothbrush with a vibration frequency of 16,000 times per minute. High-frequency vibration can easily remove gems and impurities around the teeth. It can cover all parts around the teeth surface. The sound energy generated by the motor can be transmitted to the surface of the teeth and gums through the toothbrush hairs. The adhesion of tartar and fine calculus to teeth destroys the proliferation of bacteria in the gingival pockets and tooth surfaces. It also inhibits gum inflammation and gum bleeding, and prevents gum atrophy.

### FEATURES

1. 5 modes: cleaning, sensitive, massage, whitening, gum care
2. After booting, press the button to switch the mode, long press "On/Off" button to shut down
3. Reminder function every 30S when the machine is working, automatic shutdown in 2 minutes
4. When the charging indicator flashes, it reminds that the power is insufficient and needs to be charged
5. Electric toothbrush should be keep clean and dry at charging position when charging
6. Charging method: DC straight, 5V/ 1A, 3-5 hours full, can be used for approx. 120 mins

### PRECAUTIONS

1. Before using the product, please read the user manual carefully and retain the manual for future reference. This product is suitable for children over 6 years old. The following basic safety precautions should be followed when using it. Children under 3 years old are forbidden to use. Children 4 to 6 years old need parental guidance.
2. It is recommended to replace the toothbrush head every 3 months
3. If you have had oral or gum surgery in 2 months, please consult your doctor before using this product
4. If your gums bleed a lot after using this product or you will still bleed after 1 week of use, please consult your dentist
5. This product is only designed to clean teeth, gums and tongue, please do not use it for any other purpose. If you feel any discomfort or pain, please stop using this product and contact your doctor
6. Please stop using toothbrush heads that are deformed or have cowed bristles. Please replace the toothbrush heads every 3months. If there are signs of wear, replace them in advance.
7. Do not use other toothbrush heads not recommended by the manufacturer
8. Children should be instructed when brushing their teeth until they develop good brushing habits
9. Water at the charging interface should be cleaned when charging


### SAFETY INSTRUCTIONS

Please note the following during use:

1. Remove the toothbrush head for cleaning once you brushed your teeth
2. Keep the toothbrush charging port clean and free of water when charging
3. This product is suitable for people over 6 years old
4. Product waterproof level: IPX7
5. Do not allow children under 3 years old to play with this product
6. Do not use water or acidic liquid higher than 40 to dean the charging base, the power cord and plug of the host, and do not immerse the plug in water
7. If your toothpaste contains peroxide, sodium bicarbonate or other bicarbonate (a common ingredient in whitening toothpaste), be sure to thoroughly wash the toothbrush head with soapy water after each use. It can also prevent possible plastic rupture.

### TECHNICAL SPECIFICATIONS

Name and model no.	: ecomed Sonic Toothbrush ETB 100
Input voltage	: DC 5V $\bar{\text{---}}$ 1A
Battery capacity	: 500mAh, 3.7V
Power	: 2W
Stand-by time	: 4-6 months
Charging time	: 3 hours
Operation time	: 120 minutes
Pulse frequency	: 16,000/min
Water-proof level	: IPX7
Accessories	: Toothbrush heads x 2, USB cable x 1
Dimension	: Approx. 27 x 27 x 165 mm
Weight	: Approx. 76.5 g (without toothbrush head), approx. 80 g (with toothbrush head)
Article no.	: 48131
EAN no.	: 40155 88 48131 4

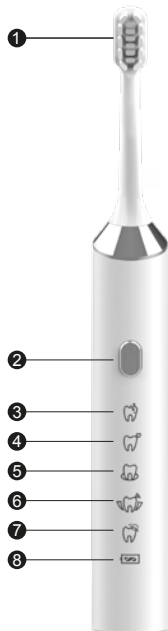- 1 Detachable toothbrush head
- 2 "On/Off" button
- 3 Clean mode
- 4 Sensitive mode
- 5 Massage mode
- 6 Whitening mode
- 7 Gum care mode
- 8 Charging indicator

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